


Evolution[®]

REVISION KNEE SYSTEM

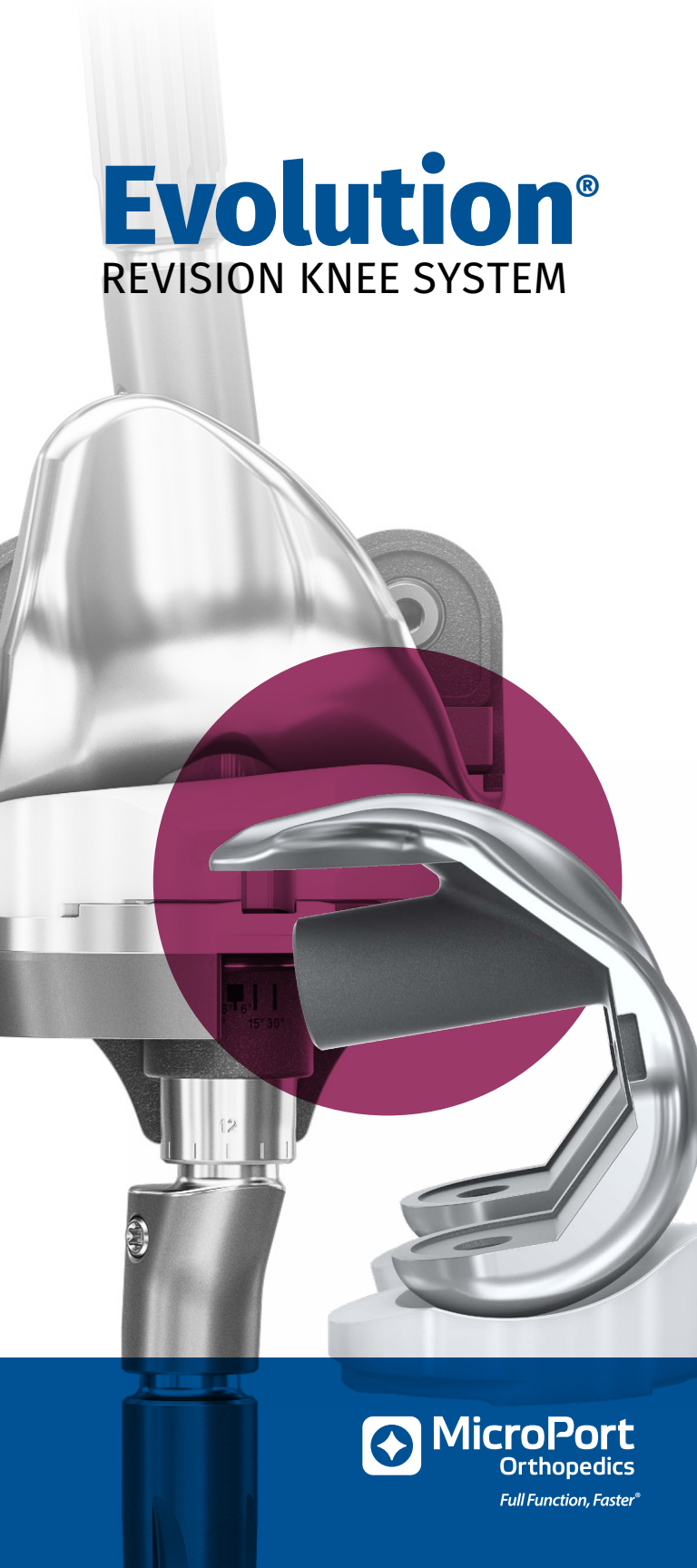

MicroPort
Orthopedics

Full Function, Faster[®]

Choose the revision option that works for your patient

FEMORAL OFFSET ADAPTERS

- 4 and 8mm offset
- 25mm length

FEMORAL AUGMENTS

- Distal and posterior
- 4, 8, and 12mm thicknesses

CCK and CS TIBIAL INSERT

- 10, 12, 14, 17, 20, 22, and 24mm thicknesses

MODULAR KEEL

- Medium and large sizes

STEM EXTENSION ADAPTERS

- 17mm diameter in 25 and 50mm lengths

CANAL-FILLING STEMS

- Diameters of 10-24mm in 1mm increments
- Lengths of 100 and 150mm

CCK and CS FEMORAL IMPLANT

- Sizes 3-8
- 35mm trunnion height

REVISION TIBIAL BASE IMPLANT

- Sizes 3-8+
- 25mm trunnion height

TIBIAL AUGMENTS

- Medial and lateral
- 5, 10, and 15mm thicknesses

CEMENTED STEMS

- 17mm diameter in 25 and 50mm lengths
- 10, 12, 14, 16, and 18mm diameters in 75mm length

The Evolution[®] Revision Knee System

Built on a 20 year, clinically-established heritage of patient satisfaction and survivorship. It is the only revision knee system that maintains the proven kinematic benefits of a medial-pivot design, while also offering surgeons intraoperative flexibility to meet all patient needs.



High medial conformity and a constant femoral radius in both planes, creates ball-in-socket articulation to maximize stability throughout range of motion




The femoral component geometry and medial-pivot insert are optimized to promote enhanced quad efficiency and all deeper flexion

Eliminate
the need
to resect
a box

The Evolution® CCK and CS inserts offer a more natural range of lateral mobility


CCK INSERT

- The lateral side of the insert has a toroid path that allows $\pm 4^\circ$ of axial rotation before cam engages with the post
- ± 2 degrees of varus valgus constraint

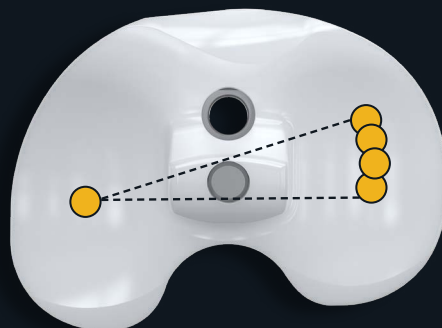

CS INSERT

- The lateral side of the insert has a toroid path that allows $\pm 15^\circ$ of axial rotation

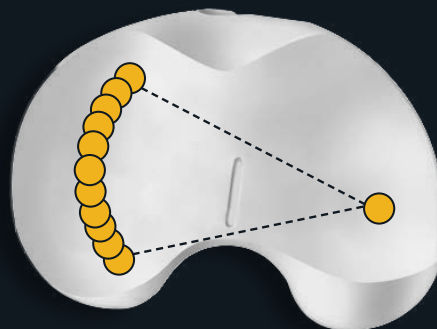



Full Function, Faster[®]



MicroPort Orthopedics Inc.
5677 Airline Road
Arlington, TN USA 38002
866 872 0211

microportortho.com

The CE-Marking of Conformity is applied per catalog number and appears on the outer package label, if applicable.

Trademarks and Registered marks of MicroPort Orthopedics Inc.
© 2019 MicroPort Orthopedics Inc. All Rights Reserved.
016904 FEB2019